

COXO[®]

Professional Dental Manufacturer

SINCE 2003

DB686 Honor

LED Curing Lights



3 Working Modes

Meeting various clinical needs

Curing Mode




Light power:
1300mw/cm

Orthodontic Mode



Light power:
>3000mw/cm²

Carise Mode

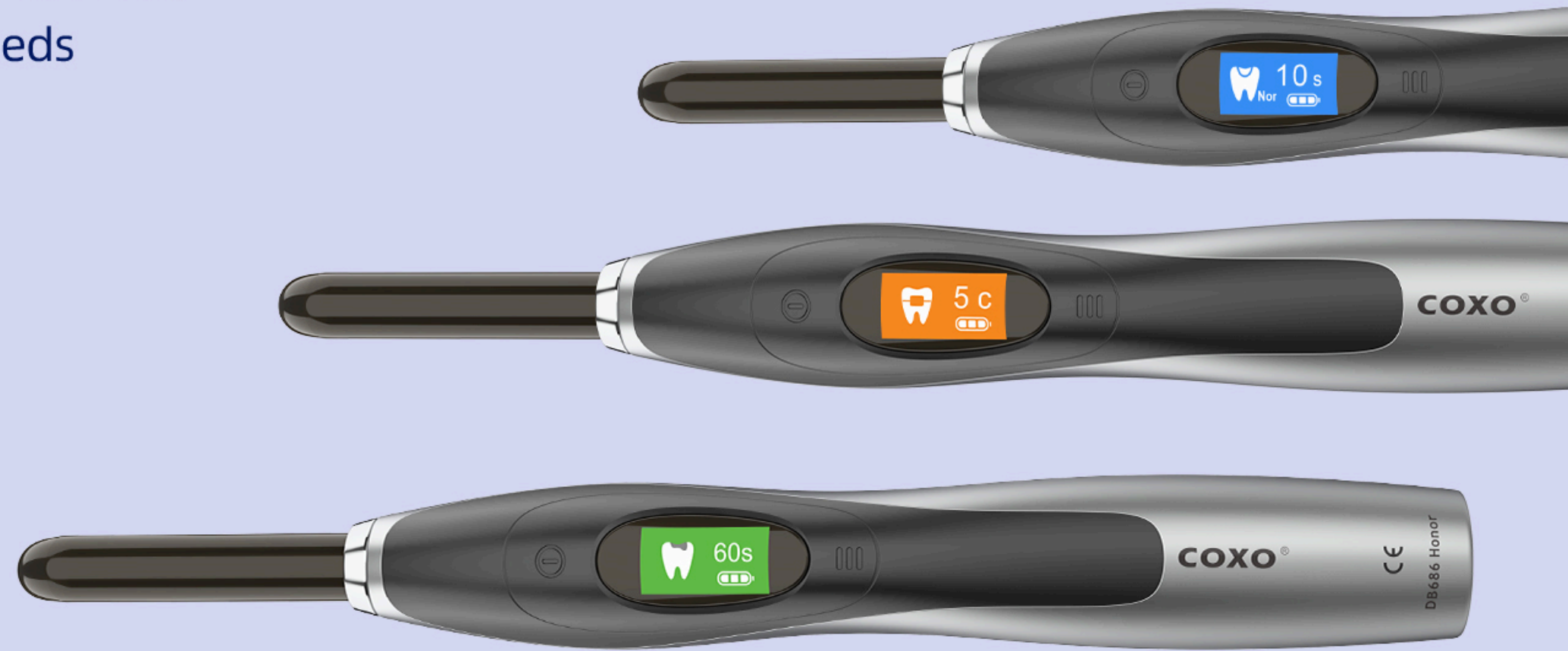


Special features



Low Temperature Curing

It can not only ensure the current effect of The resin but also received a lower temperature Rise, making patients feel more comfortable

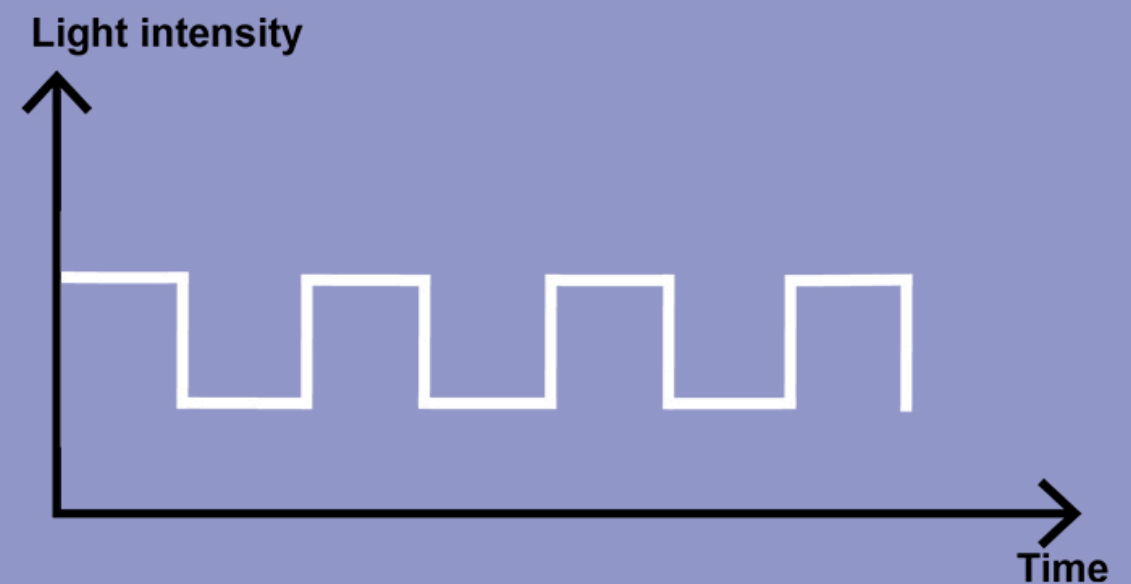

Auxiliary Focus

Spot tracking and cues can be automated to ensure our curing results are maximized



Not fully aligned, the output power will be compensated



If it's completely misaligned, the LED will not work
And it will vibrate and prompt







Wide Wavelength

385~515nm

Suitable for most of the resins in the market


Spectrum($\mu\text{W}/\text{cm}^2/\text{nm}$):0.0004
Wavelength(nm):380.0

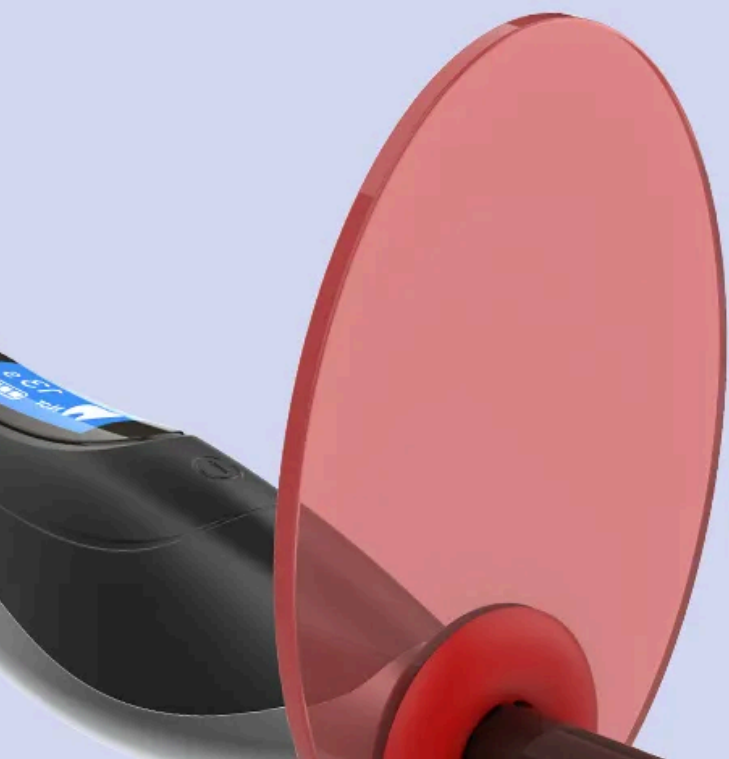

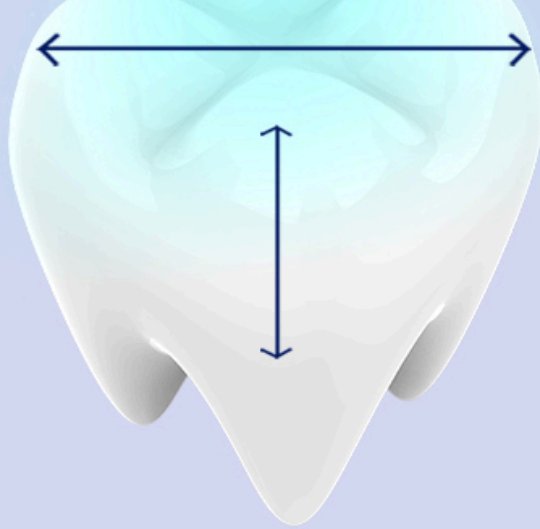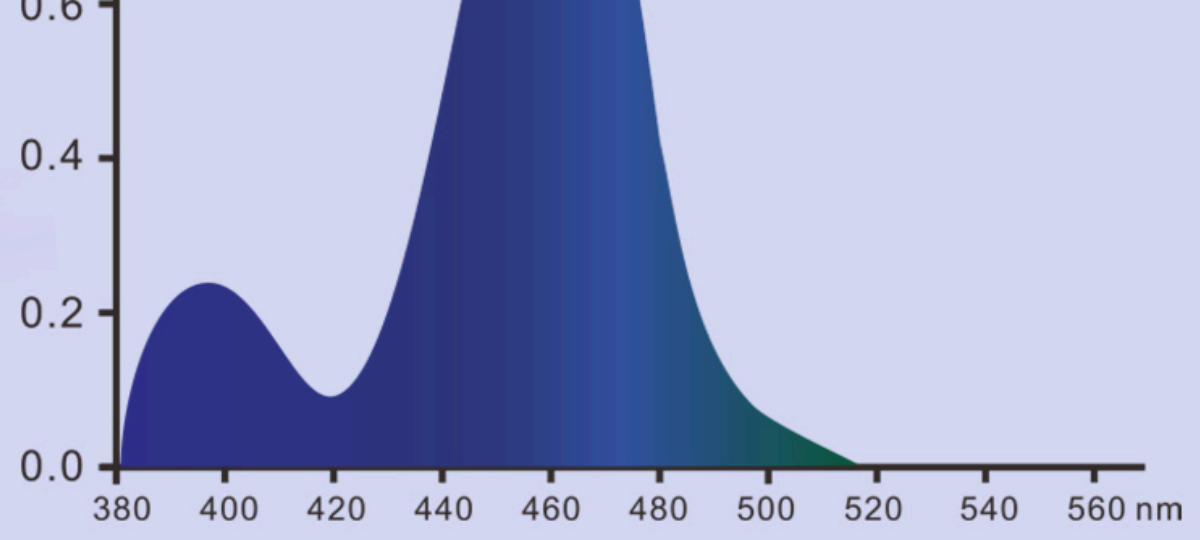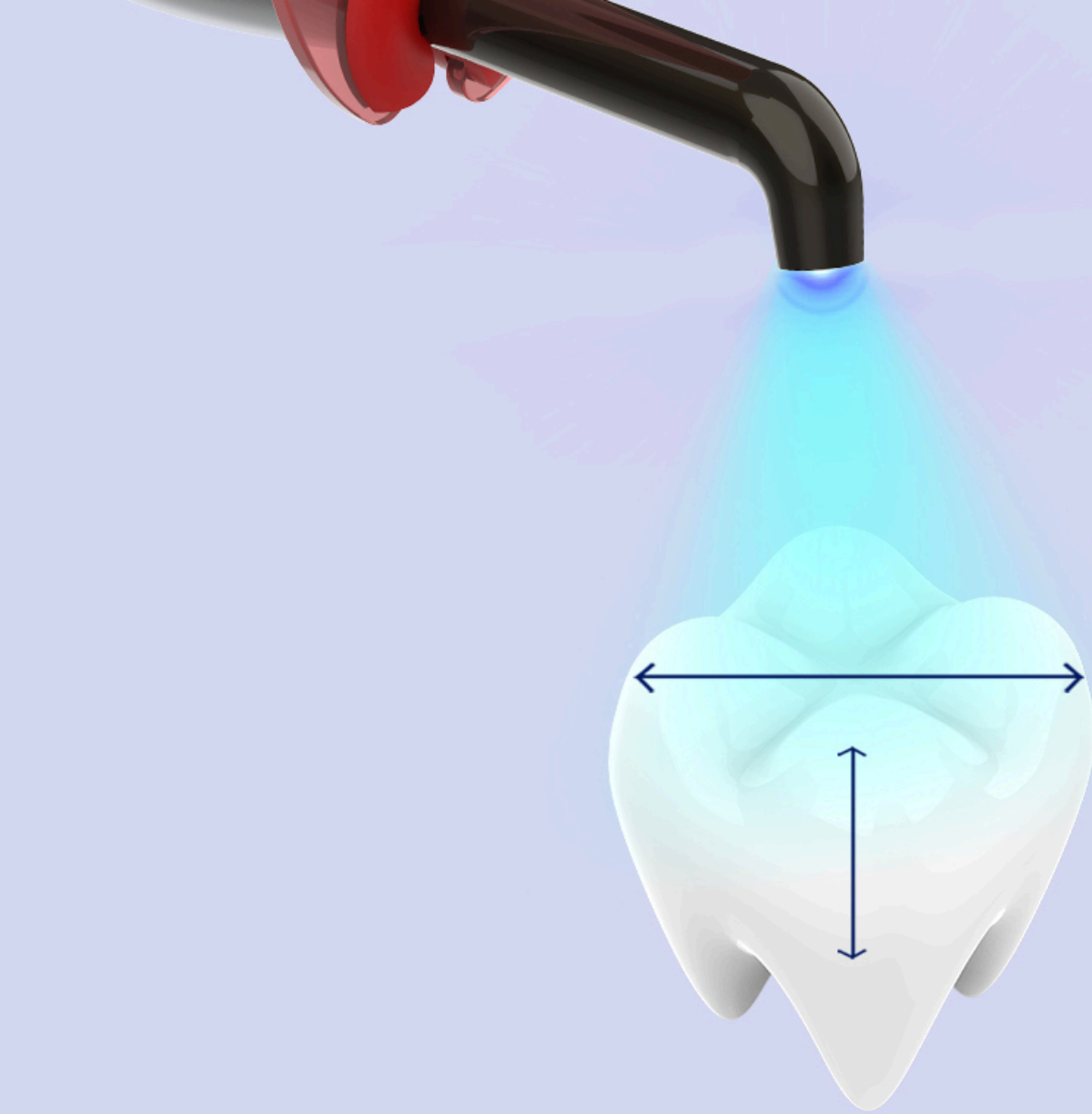

10mm large area light emission

8mm deep curing depth

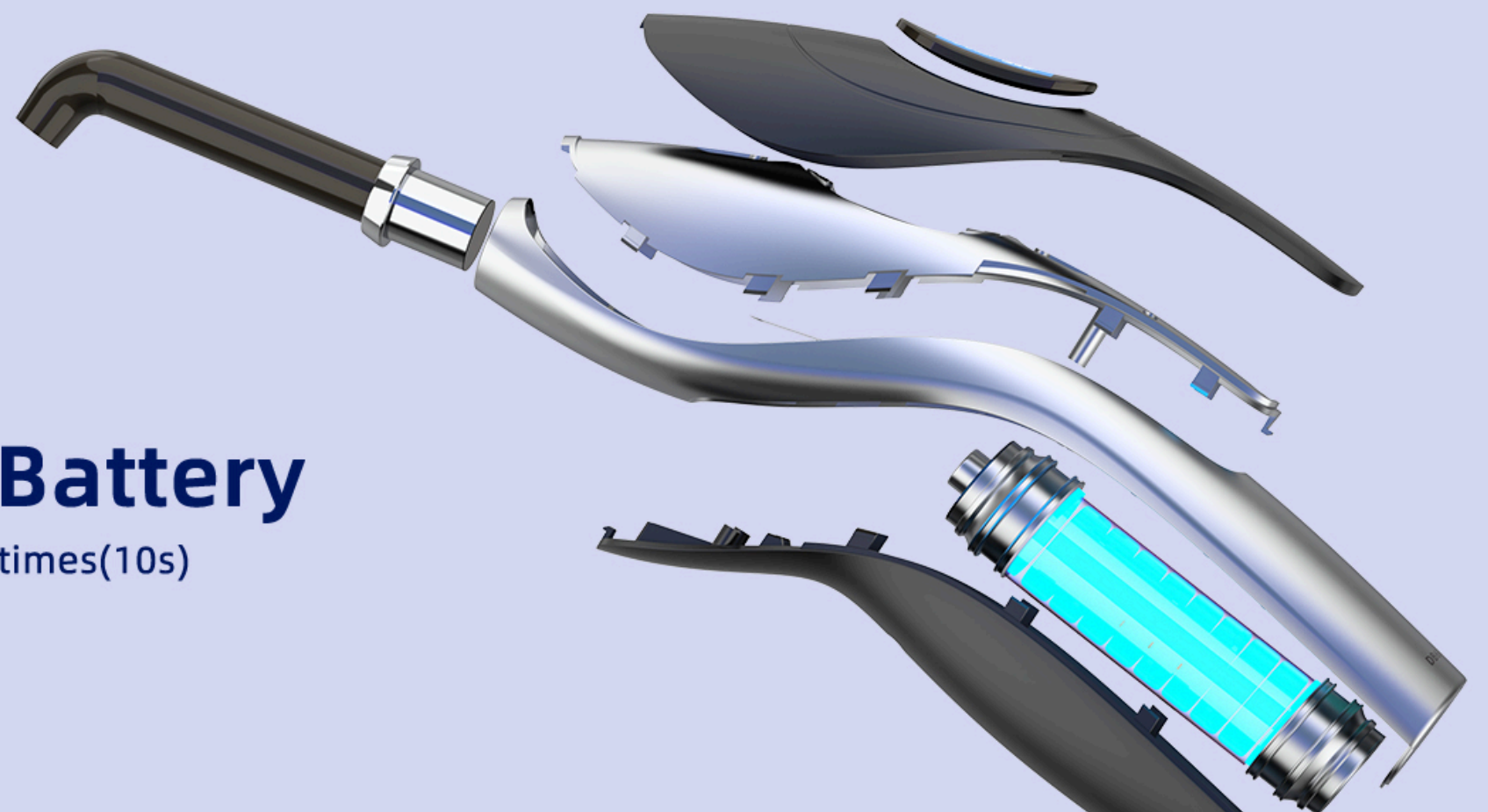


High capacity Battery

Fully charged can be used 1,000 times(10s)



Unique Grip Sensor Standby Technology

No need to repeatedly switch on/off the machine
More power saving, more humane

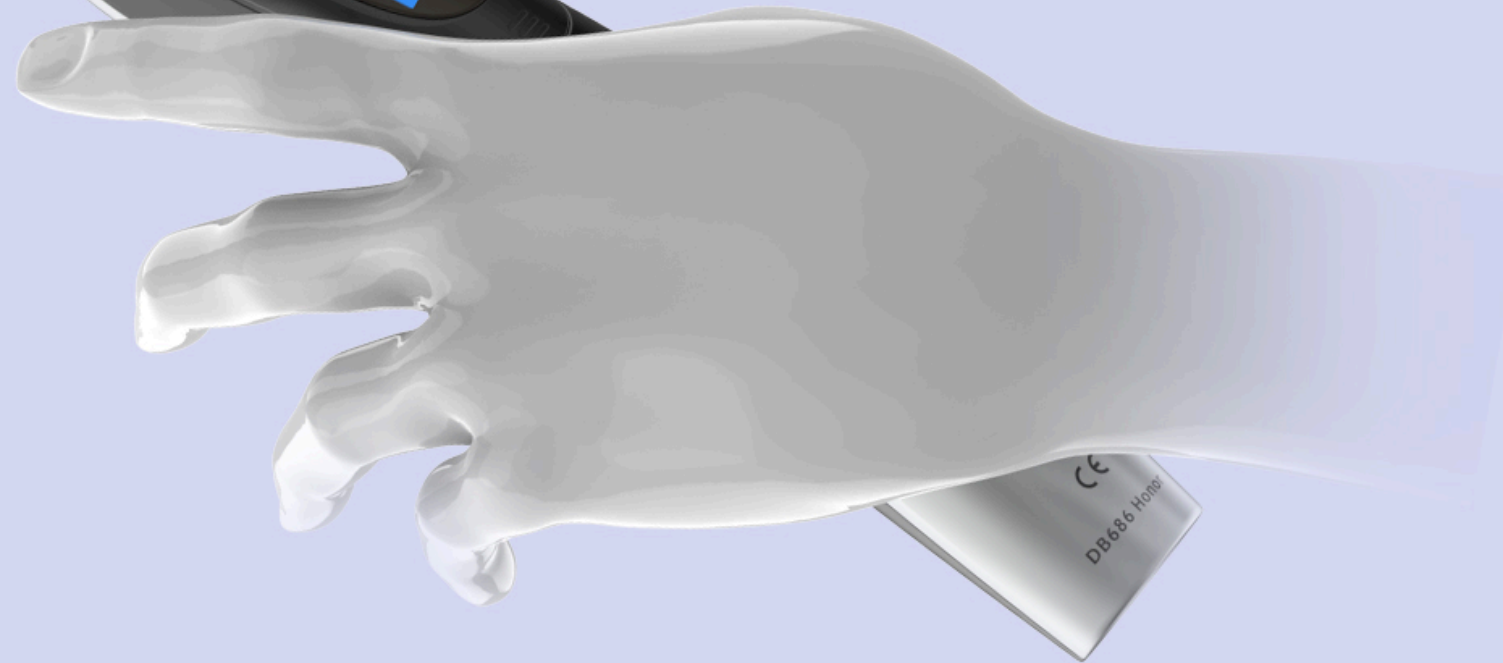

Powerful Charge Base

- Charging base with the built-in light meter
- Built-in Lithium Battery
- Non-contact recharge

